


Feature 1

Retraction can be made infinitely adjustable to allow a precise retraction.



Feature 2

Our Rochardson retractor is easy to mount on all operation tables and adjustable also for wider tables.

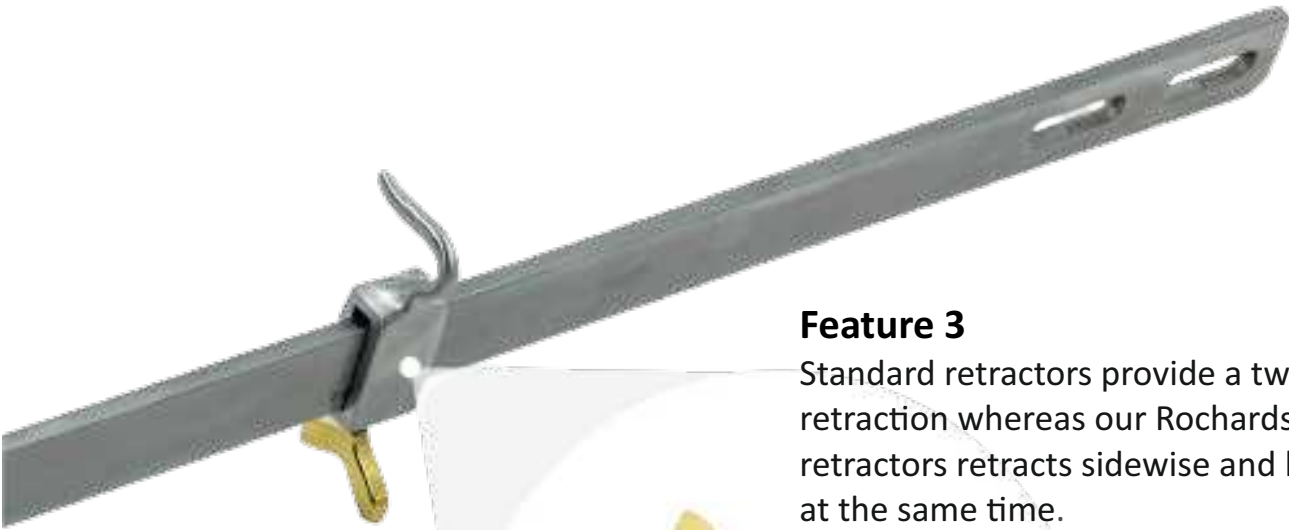

Feature 3

Standard retractors provide a two dimensional retraction whereas our Rochardson retractor retracts sidewise and lifts the thorax at the same time.



Feature 4

Blades can be adapted to each indication and patient individually, supported by a wide variety of blades.





Feature 5

The three-dimensional retraction gives excellent access even to the deeper abdomen as required for pancreas surgeries and others.

Feature 6

The retractors are manufactured of highest quality stainless steel which provides unrivalled quality and our workmanship in manufacturing the retractor is an accumulation of years of experience and quality control.





Components of Rochardson Retractor

Rail clamp TS/RR/RC/01

They are attached to the table and help in holding the side bars in place as per the desired height.



Side bar TS/RR/SB/01

They help in providing a rigid platform to hold the center rod in place. They consist of hooks on the side to hook any blades from the side profile.



Center Bar TS/RR/CB/01

The center bar is attached to the side bars with the help of a nut and bolt system keeping the center bar locked in place to help retract the body.

Center hook TS/RR/CH/01

Center hook is inserted on the center bar to hook the blades for retraction. Any number of hooks can be used as per requirement and the hooks are locked in place with the help of a wing screw.



Balfour blade 1
TS/RR/BB/01



Double angled blades
TS/RR/MAB/01



Balfour blade 2
TS/RR/BB/02 x 2 Unit

Blades

All blades are attached to a stainless steel chain to hook the blade onto the side or center bar and there is a wide range of blades to choose from.

*Blade size can be custom made as per requirement.